


CLINICAL HISTORY
MEETS INNOVATION



SAT Inc.
RESTORING MOBILITY, IMPROVING LIVES.

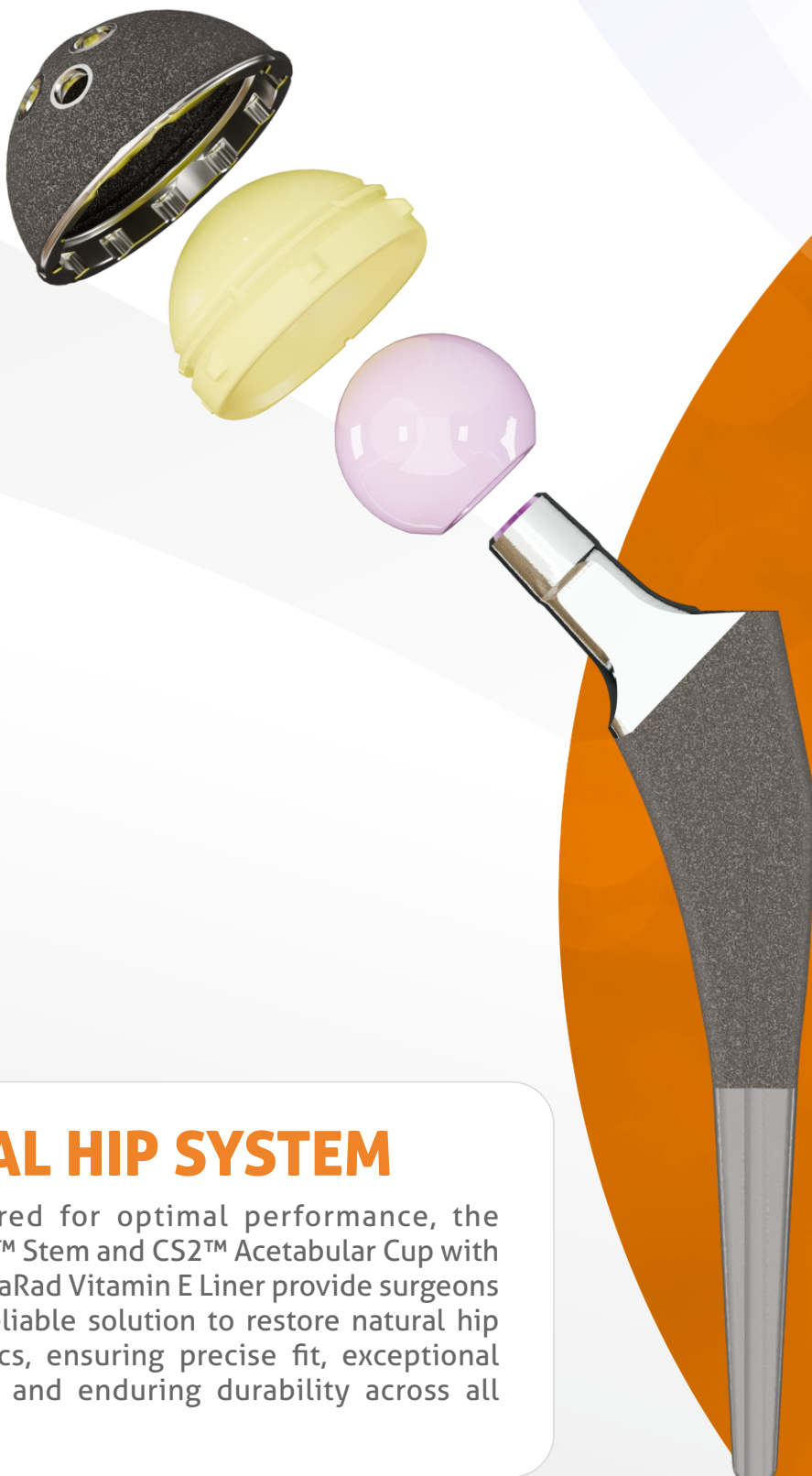

TOTAL HIP SYSTEM

Engineered for optimal performance, the Taperset™ Stem and CS2™ Acetabular Cup with a 10 MegaRad Vitamin E Liner provide surgeons with a reliable solution to restore natural hip kinematics, ensuring precise fit, exceptional stability, and enduring durability across all cases

PRECISION FIT | VITAMIN E PERFORMANCE | BIOLOGICAL INTEGRATION



TaperSet™

REDUCED DISTAL PROFILE (RDP) STEMS

The TaperSet Hip System features a reduced distal profile stem (RDP) which is designed to accommodate femurs with proximal- distal mismatch. The TaperSet RDP stem's reduced distal profile improves the proximal fit of the implant in the metaphysis.

135° NECK ANGLE

12/14 Taper Allows for restoration of joint mechanics.

DUAL TAPER WEDGE

Geometry provides stability in both mediolateral and anteroposterior planes.

POROUS PLASMA SPRAY

Proximal circumferential coating provides for biological fixation at implant-bone interface.

OFFSET OPTIONS

Standard and Lateral Offset Available

REDUCED DISTAL PROFILE (RDP)

Stems are available to address proximal distal mismatch of the femur.

Ti-6Al-4V ALLOY

Proven biocompatibility without excessive stiffness





ACETABULAR OPTIONS

CS²™ Cup with VitalitE Insert



**MULTIPLE
HOLE OPTION**



**CLUSTER
HOLE OPTION**



Vitamin E Liner

Cross-Linked Polyethylene,
Lateralized with 0°, 10°
or 20° Hoods

Enhanced Fixation

The irregular sintered bead coating enhances initial stability as well as long-term biological fixation & hemispherical design improves press-fit.


Secured Lock

Liner locking mechanism allows 12 liner positions to maximize ROM and minimize dislocations.

VitalitE CS²™ Liners


The CS2 system is available with VitalitE crosslinked poly. After 5 million cycles, Vitalite crosslinked poly has shown a 78% reduction in abrasive wear compared to cross-linked UHMWPE poly.

VITAMIN E

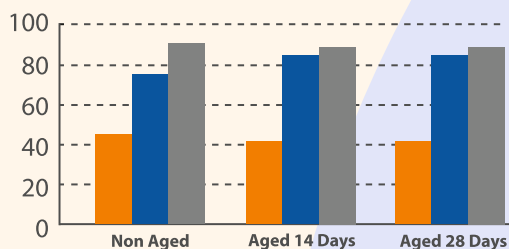

Vitalite Liners retain their long term mechanical performance

- Tensile Strength (MPa)
- Impact Strength (KJ/m²)
- Ultimate Strength (N)

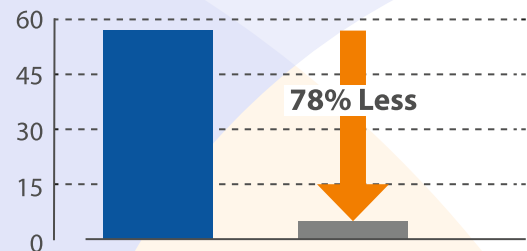


VitalitE liners demonstrated a 78% reduction in abrasive wear

- Conventional UHMWPE (mg/MC)
- VitalitE (mg/MC)



SAT Inc.
RESTORING MOBILITY, IMPROVING LIVES.

BOTH OPTIONS AVAILABLE WITH HIP SYSTEM

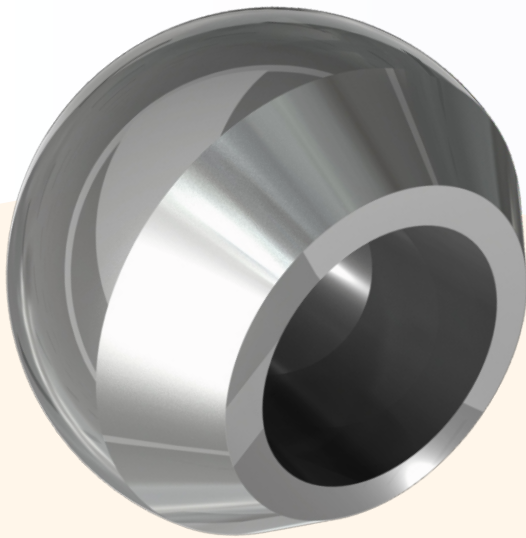

METAL HEAD



CERAMIC HEAD

For additional information, please contact your local Shalby MedTech distributor, or reach our customer service team at **+9163570 07339** or **infomedtech@shalby.org**.



SAT Inc.
RESTORING MOBILITY, IMPROVING LIVES.

USA OFFICE:
Shalby Advanced Technologies, Inc.
1115 Windfield Way,
El Dorado Hills, CA 95762.
Ph: +1916 355 7100



**SHALBY Global
Technologies Pte. Ltd.**

SINGAPORE OFFICE:
Shalby Global Technologies Pte. Ltd.
Regus One Fullerton, 1 Fullerton Rd,
#02-01, One Fullerton, Singapore 049231
Ph: +65 88707885/+60 125065062



**SHALBY
MedTech**

INDIA OFFICE:
Shalby MedTech
Opp. Karnavati Club, S.G. Highway,
Ahmedabad, Gujarat, INDIA.
Ph: +91 63570 07339