



**SAT Inc.**  
RESTORING MOBILITY, IMPROVING LIVES.

CLINICAL HISTORY  
MEETS INNOVATION



## **TAHOE UNICONDYLAR KNEE SYSTEM (TUKS)**

Simplifying **unicondylar** knee replacement surgeries with TUKS ergonomic implant and innovative instrumentation that **reduces surgery time, ensures precision and enhances clinical outcomes.**




## SIMPLIFIED & SECURE

The femoral component is a highly polished articular surface incorporating a **triple** sagittal radius and constant coronal radius for a secure **natural fit** and uniform loading conditions over a large area regardless of the flexion angle.

**R1** The sagittal and coronal radii form a true spherical surface from 10 degrees' flexion to 120 degree' flexion.

**R2** The sagittal radius becomes larger distally to minimize contact stresses in full extension.

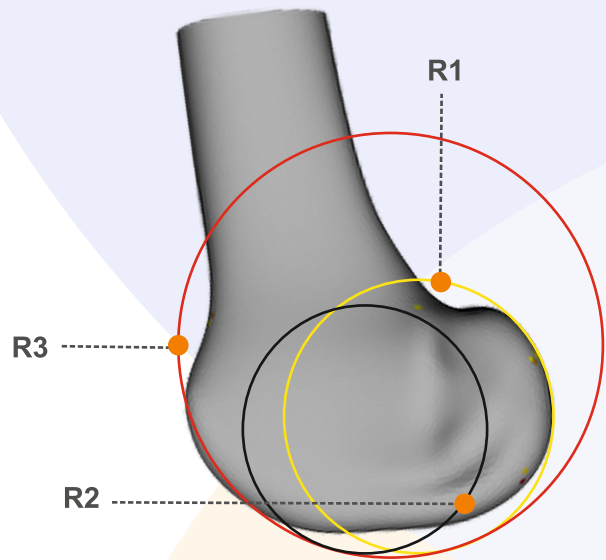



**R3** The anterior radius pulls the articular surface inward to prevent anterior impingement.

## REPRODUCIBLE

TUKS helps surgeons with efficiency and reproducibility demonstrating decreased OR times when compared to traditional instruments.

The milling technique ensures the optimum amount of bone removal leading to ideal bone resection.

TUKS instrumentation design allows a simplified surgical approach, the flexion and extension gap can now be measured at the same time.

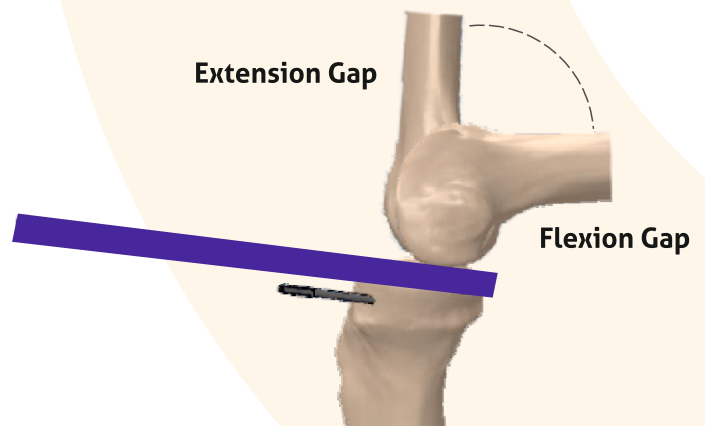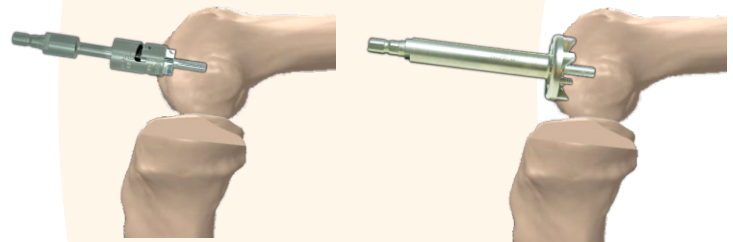



## CLINICAL

Fixed-bearing round-on-flat articulation with anatomically shaped tibial and femoral components to deliver optimal knee stability and high flexion range-of-motion

12.5% lesser bone resection and blood loss, and more bone preservation\*

VitalitE-VitalitE inserts are more resistant to oxidation than conventional cross-linked UHMWPE reducing wear by up to 80%\*\*

**TRIPLE RADIUS**  
Enables a more natural fit

**CURVED FIN**  
Enhances implant and reduces bony impact

**DUAL FEMORAL PEGS**  
Enable easy preparation & optimize fixation




**CONSTANT CORONAL RESULT**  
Allows for similar contract conditions over ROM and varus/valgus tilt


### LOCKING MECHANISM

Symmetric A/P locking and peripheral stabilizing tabs prevent pushout & micromotion

**VITALITE INSERTS**  
Vitamin E blended Inserts

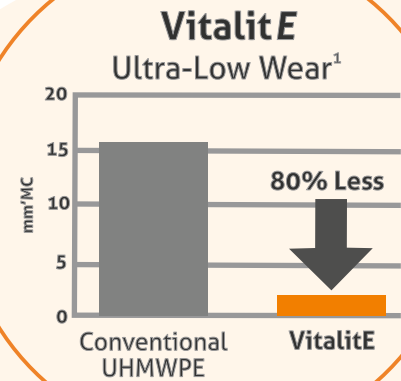

**ANATOMICAL SHAPE**  
Enables a more natural fit

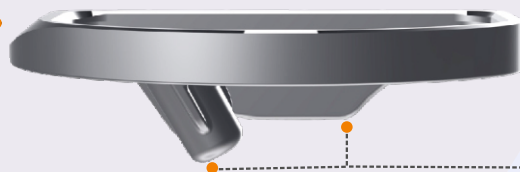

**ANTERIOR BEVEL**  
Reduces soft tissue impact

**PEG & FIN**  
Optimized alignment provides A/P & M/L stability

\*, \*\* Independent testing supported, available upon request



## INSTRUMENTATION

The innovative TUKS instrumentation allows surgeons to deliver excellence in **TUKSplasty™** Technique with its **ease of use**. The uniqueness of the instruments facilitates a **superior** way of **unicondylar** knee replacement with precise and consistent positioning of the femoral and tibial implants.



SPACER BLOCK



FEMORAL CUTTING  
BLOCK



TIBIAL PREP  
GUIDE



**SAT Inc.**  
RESTORING MOBILITY, IMPROVING LIVES.



FEMORAL STARTER  
REAMER



TRIAL FEMORAL  
IMPLANT



FEMORAL  
SPHERICAL REAMER

98119-01Rev B.1

For additional information, please contact your local Shalby MedTech distributor, or reach our customer service team at **+9163570 07339** or **infomedtech@shalby.org**.



**SAT Inc.**  
RESTORING MOBILITY, IMPROVING LIVES.



**SHALBY Global  
Technologies Pte. Ltd.**

**USA OFFICE:**

Shalby Advanced Technologies, Inc.  
1115 Windfield Way,  
El Dorado Hills, CA 95762.  
Ph: +1916 355 7100

**SINGAPORE OFFICE:**

Shalby Global Technologies Pte. Ltd.  
Regus One Fullerton, 1 Fullerton Rd,  
#02-01, One Fullerton, Singapore 049231  
Ph: +65 88707885/+60 125065062



**SHALBY  
MedTech**

**INDIA OFFICE:**

Shalby MedTech  
Opp. Karnavati Club, S.G. Highway,  
Ahmedabad, Gujarat, INDIA.  
Ph: +91 63570 07339