



SHALBY
MedTech


CLINICAL HISTORY
MEETS INNOVATION



REVSET™

Engineered to address challenges in
Revision Hip Arthroplasty

STABLE | ADAPTABLE | REPRODUCIBLE

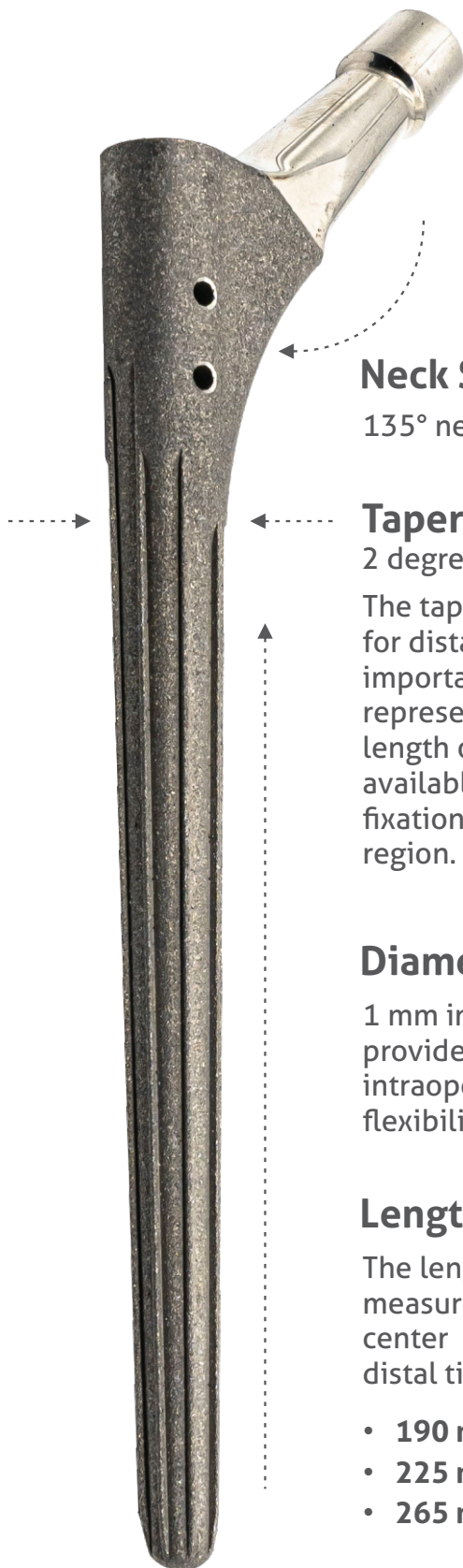

Neck Shaft Angle

135° neck shaft angle

Taper Body Stem

2 degree

The taper dimension for distally fit stem is important because it represents maximum length of implant available for distal fixation in the tapered region.

Diameter

1 mm increment
provide
intraoperative
flexibility

Length

The length is measured from center of rotation to distal tip

- 190 mm
- 225 mm
- 265 mm

3 LENGTHS



190 mm

9 size diameters
(1mm increment)
Diameters from 12 to 20mm

225 mm

9 size diameters
(1mm increment)
Diameters from 14 to 22mm

265 mm


10 size diameters
(1mm increment)
Diameters from 14 to 23mm

DESIGN FEATURES

- It is an integrated tapered stem made of Grit-blasted forged titanium alloy
- 2-degree taper of the stem body prevents subsidence and allows the anteversion angle adjustment
- Eight longitudinal ridges are designed to improve the initial stability and anti-rotational effect
- Distal tip helps with easier stem insertion
- Mirror polished neck minimises soft tissue damage & linear wear
- Grit-blasted rough surface provides potential long term stability through bone fixation
- 2 holes on proximal end (suture holes) with smooth & polished edges

REVSET's simple yet versatile instrumentation enables easy intraoperative adjustments such as additional reaming or altering stem size when necessary. It accommodates surgeon preference and intraoperative challenges.

The Swedish Hip Arthroplasty Register recommends the revision stem as the first choice.

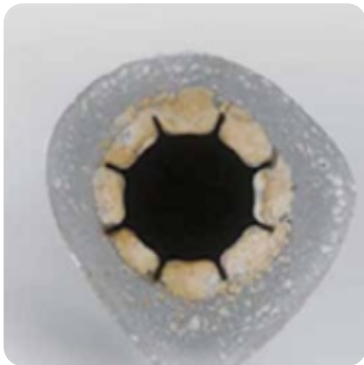

STABLE

2 degree tapered, conical implant body and eight longitudinal ribs provide rotational stability.

ADAPTABLE

REVSET offers an extensive selection of lengths and diameters, catering to varying patient anatomies and various clinical situations like artificial joint revision surgery, proximal femur osteosclerosis, bone loss, periprosthetic fractures, proximal femur deformity, difficult for proximal fixation, proximal femur fracture.



REPRODUCIBLE

Instruments are designed to maximise surgical efficiency & reproducibility of the implant position during the procedure. REVSET can be placed in any anteversion orientation by the surgeon.



SHALBY MedTech

For additional information, please contact your local Shalby MedTech distributor, or reach our customer service team at **+9163570 07339** or **infomedtech@shalby.org**.



SHALBY MedTech

Shalby MedTech

Opp. Karnavati Club, S.G. Highway, Ahmedabad,
Gujarat, INDIA.

Ph: +91 63570 07339